


# NIPPY®

## NIPPY® S+

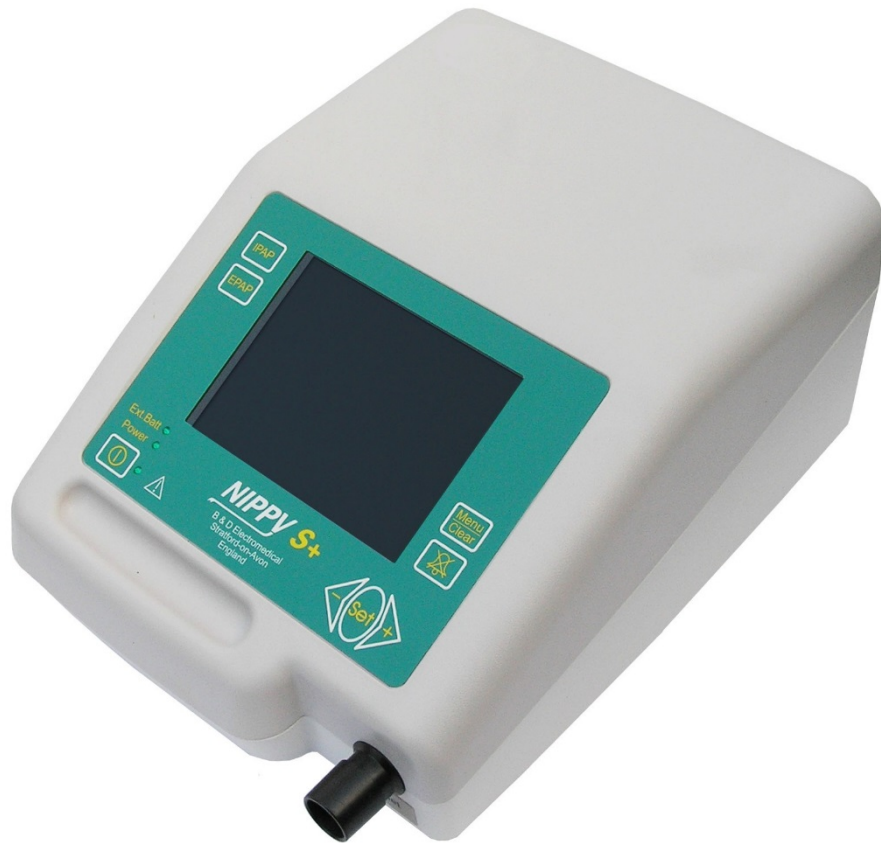

Pressure Support Ventilation

*Specialists in Ventilation & Airway Clearance*

## Pressure Support Ventilation

The NIPPY S+ is a pressure controlled positive pressure ventilator employing state of the art electronics, to give consistent and accurate support for patients requiring non-invasive ventilation via a nasal mask or nasal pillows.

## Intuitive

The large colour LCD screen gives unrivalled clarity when setting up or reading stored data. Intuitive software allows for easy setting up, with on-screen help available for each function.

Settings and airway measurements are displayed on-screen with more advanced features accessed via an easy to follow menu.

Compliance data may be viewed on-screen. Full event log may be downloaded to a PC.

The NIPPY S+ is equipped with a universal mains power inlet, a DC power inlet and an internal battery. Weighing only 3.6kg, world-wide travel has never been easier.

## Features

- Pressure Support Ventilation
- Intuitive colour LCD display
- Lockable settings
- Alarms with on-screen messages
- On-screen compliance data.
- Universal mains input and internal battery operates anywhere in the world.
- Adjustable fast response flow triggers.
- Can be connected to Nurse call systems using the optional nurse call adapter.

## Technical Specifications

Supply Voltage - 100 - 240 Vac

Dimensions (LxWxH)(cm) 30 x 22 x 13

Weight - 3.6 kg (4.5kg with battery)

Pressure Range IPAP - 3 – 38 cm H<sub>2</sub>O

Pressure Range EPAP - 3 – 20 cm H<sub>2</sub>O

Fixed Back-up Rate - 12BPM

Ti - 1.5 sec

Internal battery - 24Vdc 6.8 Ahr

Running time – 4 -12 hours

Optional External battery - 24Vdc 6.8 Ahr

Running time - 4 -12 hours

Alarms - Mask off, Power Failure